

PCTWORLD INTELLECTUAL PROPERTY ORGANIZATION
International Bureau

INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification ⁶: C12Q 1/18, 1/04	A1	(11) International Publication Number: WO 96/23901 (43) International Publication Date: 8 August 1996 (08.08.96)
(21) International Application Number: PCT/EP96/00488 (22) International Filing Date: 1 February 1996 (01.02.96) (30) Priority Data: 95200243.4 1 February 1995 (01.02.95) EP (34) Countries for which the regional or international application was filed: AT et al. (71) Applicant (for all designated States except US): GIST-BROCADES B.V. [NL/NL]; Wateringseweg 1, P.O. Box 1, NL-2600 MA Delft (NL). (72) Inventors; and (75) Inventors/Applicants (for US only): LANGEVELD, Pieter, Cornelis [NL/NL]; Venezuelastraat 55, NL-2622 BP Delft (NL). BEUKERS, Robert [NL/NL]; Sijtwinde 159, NL-2631 GZ Nootdorp (NL). BOMMELE, Michiel, Wilhelmus, Cristianus [NL/NL]; Angolastraat 24, NL-2622 EB Delft (NL). STARK, Jacobus [NL/NL]; Joost van Ospad 31, NL-3069 ZM Rotterdam (NL). (74) Agents: MATULEWICZ, Emil, Rudolf, Antonius et al.; Gist-Brocades N.V., Patents and Trademarks Dept., Wateringseweg 1, P.O. Box 1, NL-2600 MA Delft (NL).		(81) Designated States: AL, AM, AT, AU, AZ, BB, BG, BR, BY, CA, CH, CN, CZ, DE, DK, EE, ES, FI, GB, GE, HU, IS, JP, KE, KG, KP, KR, KZ, LK, LR, LS, LT, LU, LV, MD, MG, MK, MN, MW, MX, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, TJ, TM, TR, TT, UA, UG, US, UZ, VN, ARIPO patent (KE, LS, MW, SD, SZ, UG), Eurasian patent (AZ, BY, KG, KZ, RU, TJ, TM), European patent (AT, BE, CH, DE, DK, ES, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, ML, MR, NE, SN, TD, TG). Published <i>With international search report.</i> <i>Before the expiration of the time limit for amending the claims and to be republished in the event of the receipt of amendments.</i>
(54) Title: A RAPID MICROBIOLOGICAL TEST FOR THE DETECTION OF ANTIBACTERIAL COMPOUNDS (57) Abstract <p>The present invention relates to a microbiological test method for detecting antibacterial compounds in liquid samples such as milk, meat juice, serum and urine. The method employs an acid-base or redox indicator to detect growth of a test organism and is characterised in that a thermophilic strain of <i>Bacillus</i> or <i>Streptococcus</i>, e.g. <i>Bacillus stearothermophilus</i> var. <i>calidolactis</i> C953, is employed at a concentration of above 10^7 CFU/ml, preferably between 2×10^7 and 3×10^8 CFU/ml.</p>		

BEST AVAILABLE COPY

FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AM	Armenia	GB	United Kingdom	MW	Malawi
AT	Austria	GE	Georgia	MX	Mexico
AU	Australia	GN	Guinea	NE	Niger
BB	Barbados	GR	Greece	NL	Netherlands
BE	Belgium	HU	Hungary	NO	Norway
BF	Burkina Faso	IE	Ireland	NZ	New Zealand
BG	Bulgaria	IT	Italy	PL	Poland
BJ	Benin	JP	Japan	PT	Portugal
BR	Brazil	KE	Kenya	RO	Romania
BY	Belarus	KG	Kyrgyzstan	RU	Russian Federation
CA	Canada	KP	Democratic People's Republic of Korea	SD	Sudan
CF	Central African Republic	KR	Republic of Korea	SE	Sweden
CG	Congo	KZ	Kazakhstan	SG	Singapore
CH	Switzerland	LI	Liechtenstein	SI	Slovenia
CI	Côte d'Ivoire	LK	Sri Lanka	SK	Slovakia
CM	Cameroon	LR	Liberia	SN	Senegal
CN	China	LT	Lithuania	SZ	Swaziland
CS	Czechoslovakia	LU	Luxembourg	TD	Chad
CZ	Czech Republic	LV	Latvia	TG	Togo
DE	Germany	MC	Monaco	TJ	Tajikistan
DK	Denmark	MD	Republic of Moldova	TT	Trinidad and Tobago
EE	Estonia	MG	Madagascar	UA	Ukraine
ES	Spain	ML	Mali	UG	Uganda
FI	Finland	MN	Mongolia	US	United States of America
FR	France	MR	Mauritania	UZ	Uzbekistan
GA	Gabon			VN	Viet Nam

A RAPID MICROBIOLOGICAL TEST FOR THE DETECTION OF ANTIBACTERIAL COMPOUNDS

FIELD OF THE INVENTION

This invention relates to a novel microbiological test method for the rapid determination of the presence or absence of antibacterial compounds in liquids such as milk, meat juice, serum and urine. The invention also relates to a unit for the detection of residues of antibacterial compounds in liquid samples and the use of the unit.

DESCRIPTION OF THE PRIOR ART

Microbiological test methods for the determination of antibacterial compounds, particularly residues of antibiotics and chemotherapeutics such as sulphas, in liquid samples such as milk, meat juice, serum and urine have been known for a long time. Examples of such tests have been described in GB-A-1467439, EP 0005891, DE 3613794, CA 2056581 and EP 0285792. These descriptions all deal with ready to use tests that make use of a test organism and will give a result generally within 2½ to 4½, hours by the change of colour of an acid-base or redox indicator added to the test system. The principle is that when an antibacterial compound is present in the sample in a concentration sufficient to inhibit the growth of the test organism the colour of the indicator will stay the same, while when no inhibition occurs the growth of the test organism is accompanied by the formation of acid or reduced metabolites that will change the colour of the indicator.

The known test units mentioned above may include an agar medium inoculated with a suitable test organism, preferably a strain of *Bacillus* or

- 2 -

Streptococcus and an acid-base indicator or a redox indicator. The suitable test organism and the indicator are introduced into an optionally buffered agar solution, optionally nutrients are added to the solution and optionally substances to change the sensitivity to certain antimicrobial compounds in a positive or a negative way are added to the solution. Finally, the agar solution is allowed to solidify to form the agar medium in such a way that the test organisms stay alive but cannot multiply because of lack of nutrients and/or low temperature. Optionally compounds may be added to the agar medium, e.g. the nutrients or the indicator(s) as a separate source, e.g. as a tablet or a paper disc.

A preferred species of *Streptococcus* for use in conventional microbiological tests for antibacterial compounds is *Streptococcus thermophilus*. A preferred species of *Bacillus* for use in such tests is *Bacillus stearothermophilus*. This species may be employed for this purpose as spores.

In most documents concerning microbiological tests for antibacterial compounds, a very broad concentration of the test organism of 10^5 to 10^8 - 10^9 colony forming units (CFU) per ml is mentioned. Examples of such documents are GB-A-1467439, EP 000581, EP 0611001, DE 3429823, US 3941658, US 4946777 and EP 0285792. However, none of these documents has an example specifying a concentration of test organism higher than 10^7 CFU/ml. This reflects the fact that in practice microbiological tests for antibacterial compounds have previously routinely been carried out with a test organism concentration below 10^7 CFU/ml. In CA 2056581 (page 11, line 27), it is specifically stated that the concentration of test organism - in this particular case the spores of *Bacillus stearothermophilus* - may vary, but may never exceed a concentration of 10^7 CFU per ml agar medium.

Thus, up till now, it has been commonly accepted that a test organism concentration of 10^7 CFU per ml is the maximum concentration that could be used in the type of microbiological tests described in this document.

Up till now, it has been commonly thought that use of a higher microorganism

concentration in such a test will result in loss of sensitivity for relevant antibacterial compounds.

Of course, a suitable rapid test for antibacterial compounds should among others fulfil the following requirement:

- 5 -a high sensitivity for a wide range of antibacterial compounds used in practice.

A test with a bacterium concentration of 10^6 to 10^7 CFU per ml generally gives a result between 2½ to 4½ hours or even after a longer time. Although a test duration of minimal 1½ hours sometimes is mentioned (for
10 example EP 000581, GB-A-1467439, US 3941658, US 4946777) such a test has never been demonstrated in an example or anywhere else.

Up till now, it was not considered possible to carry out a microbiological test for antibacterial compounds with a test duration less than 2½ hours by using a concentration of test organism higher than 10^7 per ml
15 while retaining good sensitivity for a wide range of antibacterial compounds. The present specification now teaches that, contrary to expectation, such microbiological tests are feasible.

DESCRIPTION OF THE INVENTION

20

In one aspect, the present invention provides a method for detecting one or more antibacterial compounds in a liquid sample, such as milk, meat juice, serum or urine, which comprises the steps of:

- (i) contacting said sample with therein a test microorganism selected
25 from thermophilic strains of *Bacillus* and *Streptococcus*, prior to addition of said sample growth of said test microorganism being prevented;
- (ii) providing conditions for said test microorganism whereby in the absence of antibacterial compounds growth of said test microorganism can occur and can be detected; and
- 30 (iii) detecting the growth of said microorganism after a time sufficient for detection of test microorganism growth in the absence of antibacterial

compounds, but during which period detectable test microorganism growth is maintained by a known amount of an antibacterial compound in a test sample; characterised in that said test microorganism is present at above 10^7 CFU/ml.

Such a test method provides a marked reduction in test duration
5 compared with conventional microbiological tests for antibacterial compounds. While conventional microbiological tests for antibacterial compounds generally give a result between $2\frac{1}{2}$ and $4\frac{1}{2}$ hours, a method described in this document gives a result in less than $2\frac{1}{2}$ hours, e.g. 70 minutes or possibly less e.g. 60 minutes or less. This is of extreme importance to the user because the quality
10 of the sample liquid is known more quickly thus allowing an earlier delivery o. processing.

Moreover, surprisingly, by employing a microbiological test of the invention with a test organism concentration of more than 10^7 CFU per ml, a result can be obtained in less than 150 minutes without losing sensitivity. Thus,
15 it has been found possible using a microbiological test of the invention to detect in well below $2\frac{1}{2}$ hours, e.g. as little as 70 minutes, Penicillin G at 0.003 IU/ml in a milk sample (see Examples).

The test microorganism present in the test is preserved under such conditions that they are activated, even after at least one year, when contacted
20 with the test sample.

The test microorganism may be preserved in dispersed form in a solidified medium e.g. solidified agar medium.

According to another embodiment the test microorganism is preserved in concentrated form e.g. in dried form. In case the microorganism is dispersed
25 in solidified medium, the amount or presence of the microorganism is defined as CFU of said microorganism per ml of said medium (CFU/ml). In case the microorganism is present in other forms (e.g. freeze-dried), the amount or presence of the microorganism is defined as CFU's of said microorganism per ml test sample (e.g. milk) added during the method of detection (CFU/ml).

30 Antibacterial compounds which can be detected in accordance with the present invention are antibiotics and chemotherapeutics like sulpha

- 5 -

compounds. A test method of the invention may be used for the detection of antibacterial compounds in a wide range of liquid samples, for example milk, milk products, meat juices, serum, urine, eggs, honey, water and extracts of meat, eggs or honey.

5 In a further aspect, the invention provides a test unit for carrying out a method of the invention comprising a container holding a solidified agar medium having dispersed therein a test organism selected from thermophilic strains of *Bacillus* and *Streptococcus* at above 10^7 CFU/ml, said agar medium optionally additionally containing nutrients for growth of said test organism
10 and/or one or more components for detection of said growth. Examples of containers useful for the purpose of the invention are tubes or blocks with tubular indentations, for example, transparent tubes, single or in a set, or combined in the form of a block of translucent material having a number of holes shaped therein, e.g. microtiter plates.

15 The agar medium in such a unit may be an optionally buffered agar medium to which is added a test organism, a nutrient source comprising sufficient nutrients to support growth of the test organism, an acid-base indicator or one or more redox indicators and optionally substances to change the sensitivity to certain antimicrobial compounds in a positive or a negative
20 way. All the test components may be incorporated into the initial agar solution used for preparation of the test unit or at least some may be added separately to the agar medium.

Thus, optionally at least part of the nutrients may be subsequently incorporated into the agar medium in the test unit from a separate applied
25 source, for example a tablet or paper disc, placed on the agar medium before carrying out the test. At least one indicator component may also be added separately to the agar medium, e.g. by inclusion in a separate nutrient source. Also a buffer and the test organism may be added as a separate source. The nutrient source and / or one or more indicator components, e.g. redox
30 indicators, may be added to the test sample.

- 6 -

When nutrients and / or indicator (s) and / or the test organism are added to the agar medium or test sample from a separate applied source, e.g. a tablet, measures are preferably taken to avoid moisture transport from the agar medium into the applied additive source. When a tablet is employed for this purpose, this may be done by coating the tablet with a moisture resistant layer, for example a wax, which coating must degrade or melt during the test procedure. A wax having a melting temperature of 35° to 55°C, preferably 40° to 45°C is suitable.

The nutrients are added to enable the multiplication of the test organism in the absence of antibacterial compounds. Suitable nutrients are for example assimilable carbon sources (e.g. lactose, glucose or dextrose) , assimilable nitrogen sources (e.g. peptone) and sources of growth factors and minerals (e.g. yeast extract).

The growth of the test microorganism can be detected by colour change of the agar medium or the test sample resulting from the presence of an acid base or redox indicator.

If the indicator used is an acid-base indicator, an indicator for a pH of about 5.5 is suitable, preferably bromocresolpurple or phenolred. Other acid-base indicators may however be used as well.

If it is desired to use a redox indicator or a combination of two or more redox indicators, suitable indicators are e.g. Brilliant Black, Methylene Blue, Toluidine Blue, Nile Blue A, 2,3,5-triphenyltetrazolium, Safranin O, Indigo Carmin, Thionin, Gallocyanine, Brilliant Crocein MOO, Acid Yellow 38, Acid Orange 51, Acid Blue 120, Basic Blue 3, Azure A, Azure B, Congo Red, 1-10 Phenanthroline, Janus Green B, Brilliant Cresyl Blue. Other redox indicators (redox mediators, redox catalysts and electron carriers) may be used as well. Preferred redox indicators are Brilliant Black, Methylene Blue, Toluidine Blue and Nile Blue A or combinations thereof, e.g. a combination of Brilliant Black and Toluidine Blue.

- 7 -

A suitable test organism for a method of the invention is *Streptococcus thermophilus*, preferably *Streptococcus thermophilus* T101 (DSM 4022, deposited on March 3, 1987).

Bacillus stearothermophilus may alternatively be employed. Particularly preferred for use in a method of the invention is *Bacillus stearothermophilus* var. *calidolactis* C953. Strain C953 of *Bacillus stearothermophilus* var. *calidolactis* was deposited with the Laboratory of Microbiology of the Technical University of Delft under the accession number LMD 74.1 in 1974 and with the Centraal Bureau voor Schimmelcultures (CBS), Baarn under the accession number CBS 760.83 in 1983 where the strain is available to the public.

Bacillus stearothermophilus var. *calidolactis* C953 and *Streptococcus thermophilus* T101 are very sensitive to antibacterial compounds, especially antibiotics such as penicillins and chemotherapeutics such as sulpha compounds. These microorganisms are fast growing and have the additional advantage that their optimal growing temperature is high (*Bacillus stearothermophilus*: between 50°C and 70°C; *Streptococcus thermophilus*: between 35°C and 45°C). This makes these microorganisms very suitable for a test according to the invention as they cannot grow at storage temperature (e.g. room temperature) and there is little possibility that organisms which are possibly present in the test liquid or which have otherwise been included in the test system would affect the result of the test.

When the test organism is a *Bacillus* strain, it is preferably incorporated into the agar medium in the form of a spore suspension which may be prepared and incorporated into the agar medium prior to solidification by known methods (see, for example, GB-A-1467439). Alternatively, as hereinbefore indicated, the test organism may be added as a separate source, for example by application of a tablet or paper disc to the agar medium.

When the test organism is a *Streptococcus* strain, the bacteria are preferably incorporated into the agar medium in the form of bacterial cells which may be prepared according to known methods (EP 0 285 792).

When the test microorganism is in concentrated form (e.g. freeze-dried form, tablet, paper disc, etc.) the concentrate is measured into a vessel which may be a conventional ampoule, a sealable test tube, a sample bottle, a hole of a block such as a microtiter plate and the like (see for example EP-0285792).

5 The concentration of test organism in the agar medium is suitably above 10^7 CFU/ml, but desirably does not exceed 10^{10} CFU/ml, e.g. a test organism concentration between 2×10^7 and 10^9 CFU/ml may be employed, preferably 2×10^7 to 3×10^8 CFU/ml. Most preferably, the test organism concentration in the is between 3×10^7 and 2×10^8 CFU/ml.

10 According to a preferred embodiment of the present application, the sensitivity of the test organism is adjustable. The sensitivity may be altered by various means, for example by adding certain substances, by changing the test conditions such as the pH or concentration of buffering substances or agar or by varying the ratio of the volumes of agar and test liquid.

15 Examples of substances that may be added to the test system to change sensitivity are nucleosides such as adenoside and antifolates such as trimethoprim, ormethoprim and tetroxoprim which improve the sensitivity of the test organism to sulpha compounds. Salts of oxalic acid or hydrofluoric acid may be added to improve the sensitivity to tetracyclines. Cysteine may be
20 added to diminish the sensitivity to penicillins.

Hereinafter the invention will be discussed in more detail for a test comprising the agar medium, however, the skilled in the art will appreciate that these teachings can also be applied for tests without an agar medium.

As hereinbefore indicated, components needed or optionally desirable in
25 the agar medium, including optionally substances to change the sensitivity to certain antibacterial compounds, may be added to an optionally buffered agar solution. The agar solution is allowed to solidify in the test units to form the agar medium in such a way that the test organism stays alive but cannot multiply because of lack of nutrients and / or low temperature. The needed
30 compounds which have not been added to the agar solution prior to

solidification have to be added from a separate source, for example a tablet or paper disc.

5 Needed components are the test organism, one or more indicators and the nutrients. Other compounds, e.g. substances to improve or diminish the sensitivity to certain antibacterial compounds and buffer solutions, are optional.

10 It is preferred to carry out the process of the present invention in such a way that the test organism does not grow or multiply in the agar medium in the presence of a known minimum amount of an antibacterial compound to be detected. Also the test units, e.g. test tubes containing the agar medium, have to be stored in such a way that the test organism, does not multiply in the agar medium. This is generally achieved by depriving the organism of nutrients until the test is carried out or by maintaining the culture (during the production process) and the test units (during storage) at a sufficiently low temperature, preferably between 4°C and 30°C, more preferably between 4°C and 15°C, or
15 both.

The test units are preferably closed air-tight during storage in which condition they may be stored for at least several months.

20 The optionally buffered agar medium containing the test organism and / or the indicator(s) and / or the nutrients and / or certain substances to alter the sensitivity of the test is allowed to solidify in upright test units e.g. tubes.

The units preferably have determined sizes. This is because of the reliability of the test. The height of the agar medium in the test units is of great importance. If an antibiotic is present in the test sample it will diffuse into the agar. The sensitivity of the test is therefore preferably ascertained with variation
25 of the height of the solidified agar in the test unit.

Preferably, the height of the agar medium in the test unit is selected in such a way that the indicator changes its colour when the antibacterial compound to be detected is present in a concentration below a certain value, but does not change its colour when the concentration is above this value.
30 Preferably, the test is arranged to be a test from which one can read whether a certain concentration of an antibacterial compound is present or not. For this,

- 10 -

the height of the test unit will preferably be high enough to contain an amount of agar medium and sample corresponding to a height of 3-30 mm, more preferably 5-15 mm. The internal cross-sectional dimension of the test units is preferably 1-20 mm, more preferably 1-14 mm.

5 The volume of the agar medium in the test unit is determined by the height of the test unit, the internal cross-sectional dimension of the test unit and the percentage of the volume of the test unit which is filled with the agar medium. The volume of the agar medium is preferably 10 μ l - 3 ml, more preferably 15 μ l - 1 ml and most preferably 20 μ l - 500 μ l.

10 A suitable liquid test sample volume is for example 0,01 - 0,5 ml. Thus, for example, a test unit of the invention formed by solidifying 0.3 ml agar solution in an upright sterile test tube of internal diameter 9 mm has been found suitable for detecting 0.003 IU/ml Penicillin G in a milk sample.

 When spores of *Bacillus stearothermophilus* var. *calidolactis* are used,
15 suitable incubation temperatures are between 55°C and 70°C, more preferably between 60°C and 66°C. Incubation can be carried out in any conventional form of incubator, e.g. a waterbath, a block heater, a stove or a dry incubator.

 When cells of *Streptococcus thermophilus* are used, suitable incubation temperatures are between 36 and 44°C, more preferably between 38 and 42°C.

20 As hereinbefore indicated, an incubation period of less than 150 minutes is feasible, e.g. 70 minutes, possibly as little as 60 minutes or even 30 minutes. Such an incubation period is markedly shorter than the incubation period of other conventional microbiological test systems for antibacterial compounds employing less than 10⁷ CFU's/ml.

25 A test of the invention is very simple to carry out, so that persons who perform the test do not have to be specially educated. The result is very simple to determine. After the incubation period, the colour of the agar medium containing the indicator shows if test organism growth did occur or not.

 All documents mentioned in this application are herein incorporated by
30 reference to the same extent as if each individual application or patent was specifically and individually indicated to be incorporated by reference.

Legend to the Figures

Fig. 1: shows the test duration as a function of the spore concentration for test units of type A employed as in Example IV.

5 Fig.2: shows the test duration as a function of the spore concentration for test units of type B employed as in Example V.

EXAMPLE I

Preparation of test tubes of type A to detect antibiotics

10

Test tubes of type A contained an agar medium prepared as follows:

A solution was made of 12 g agar, 9 g sodium chloride and 50 ml of a 0,1 M triethanolamine buffer solution (pH = 8.0) in 1000 ml distilled water. An amount of a solution of bromocresolpurple to give a final concentration of 20 mg per 1000 ml was also added. The final solution was sterilized for 20 minutes at 121°C and cooled to about 60°C.

15

To this solution different amounts of a suspension of *Bacillus stearothermophilus* var. *calidolactis* C953 spores in distilled water were added to give final concentrations between 10^6 and 10^{10} spores per ml.

20

Sterile tubes with a diameter of about 9 mm were filled with 0,3 ml of the agar solution under aseptic conditions and immediately sealed, e.g. with aluminium foil. The agar solution in the test tubes was allowed to solidify in the test tubes while the test tubes were held in an upright position.

The tubes were stored at a temperature between 5°C and 15°C.

25

EXAMPLE II

Preparation of test tubes of type B to detect antibiotics

Test tubes of type B contained an agar medium prepared as follows:

30 A solution was made of 10 g agar, 9 g sodium chloride, 2 g glucose, 2 g peptone, 2 g tryptone, 4 g yeast extract and 50 ml of a 0.1 M triethanolamine

- 12 -

buffer solution (pH = 8.0) in 1000 ml distilled water. The final solution was sterilized for 20 minutes at 121°C and cooled to about 60°C.

An amount of a solution of Brilliant Black to give a final concentration of 80 mg per 1000 ml and an amount of a solution of Toluidine Blue to give a final concentration of 3 mg per 1000 ml were also added.

To this solution different amounts of a suspension of *Bacillus stearothermophilus* var. *calidolactis* C953 spores in distilled water were added to give final concentrations between 10^6 and 10^{10} spores per ml.

Sterile tubes with a diameter of about 9 mm were filled with 0,3 ml of the agar solution under aseptic conditions and immediately sealed e.g. with an aluminium foil. The agar solution in the test tubes was allowed to solidify in the test tubes while the test tubes were held in an upright position.

The tubes were stored at a temperature between 5°C and 15°C.

15 EXAMPLE III

Preparation of nutrient tablets

A mixture was made of 100 g dextrose, 160 g Avicel PH101, 50 g tryptose, 14 g yeast extract and 15 g precinol. This mixture was sufficient to prepare about 18000 tablets with a diameter of approximately 3 mm and a thickness of approximately 1.9 mm.

EXAMPLE IV

25 Carrying out the test type A

One of each of the test tubes of type A with different concentrations of *Bacillus stearothermophilus* var. *calidolactis* spores as described in Example I were opened by removing the seals. A nutrient tablet as described in Example III was added to each of the test tubes. Then 0.1 ml of a fresh cow's milk

- 13 -

sample was added to each of the test tubes and the test tubes were immediately placed in a waterbath kept at 64°C.

Observations were made after 30 minutes to 2 hours and 30 minutes. The time at which the colour of the agar medium became yellow was determined. As is shown in Figure 1, the test duration correlated with the
5 number of spores added to the agar medium.

If at this specific time, the colour of the agar medium was blue, the sample was known to contain a detectable amount of an antibacterial compound (e.g. 0.003 I.U. per ml or more penicillin G). Figure 1 shows that a
10 test duration reduction to at least 60 minutes is realizable.

EXAMPLE V

Carrying out the test type B

15 The procedure described in example IV was followed except that test tubes of type B were used and no tablet was added to the test tube.

Figure 2 shows that a test duration reduction to at least 70 minutes is realizable.

20 EXAMPLE VI

Carrying out the test type B

An amount of penicillin G of 0.003 I.U. per ml was added to fresh cow's milk. The control sample contained no penicillin G.

25 Two of each of the test tubes of type B with concentrations of *Bacillus stearothermophilus* var. *calidolactis* spores as described in Example II were opened
by removing the seal. Then 0,1 ml of each of the two milk samples was added to each of the test tubes and the test tubes were immediately placed in a
30 waterbath kept at 64°C.

Observations were made after 30 minutes to 3 hours. The time at which the colour of the agar medium became yellow was determined. As is shown in Table 1, the test duration correlated with the amount of spores added to the agar medium without influencing the sensitivity of the test in a negative way.

- 5 Table 1 shows a decreasing test duration down to at least 70 minutes. Even at a test duration of 70 minutes, a concentration of 0.003 I.U. per ml penicillin G can be detected. Of course also higher concentrations of penicillin G can be detected.

- 15 -

Table 1 Sensitivity of test type B at different spore concentrations

test duration (minutes)	180	155	140	130	120	110	100	70
spores/ml medium	$6,5 \times 10^6$	$1,3 \times 10^7$	$2,0 \times 10^7$	$2,6 \times 10^7$	$3,3 \times 10^7$	$4,9 \times 10^7$	$6,5 \times 10^7$	$1,3 \times 10^8$
milk without antibiotic	-	-	-	-	-	-	-	-
milk with 0.003 IU/ml pen.	+	+	+	+	+	+	+	+

EXAMPLE VIIDetection of antibiotics and sulpha compounds in microtiterplates (test type C)

5 For the preparation of the test plates (test type C) the procedure described in Example I was followed except that the nutrients described in Example III (dextrose, tryptose and yeast extract) were not added separately but were added to the initial agar medium in the same concentration as described in Example I. Also together with the buffer solution an amount of trimethoprim
10 solution was added to give a final concentration of 60 µg per litre.

To the wells of sterile microtiterplates 20 µl of the agar solution was added under aseptic conditions. The plates were immediately sealed with aluminium foil.

The plates were stored at a temperature between 5°C and 15°C.

15 An amount of penicillin G of 3 ppb was added to fresh cow's milk. Also a milk sample was prepared which contains an amount of sulphadiazine of 100 ppb.

The control sample contained no penicillin G or sulphadiazine.

The microtiterplates were opened by removing the seal. Then 20 µl of
20 each of the three milk samples was added to each of the test plates (duplo) containing different concentrations of *Bacillus stearothermophilus* var. *calidolactis* spores C953. The plates were immediately placed in a waterbath kept at 64°C.

Observations were made after 30 minutes to 150 minutes. The time at
25 which the colour of the agar medium became yellow was determined. As is shown in table 2, the test duration correlated with the amount of spores added to the agar medium without influencing the sensitivity of the test in a negative way.

Even at a test duration of 90 minutes both penicillin G and
30 sulphadiazine can be detected at concentrations of respectively 3 ppb penicillin

- 17 -

G and 100 ppb sulfadiazine. Penicillin G (3 ppb) can even be detected at a test duration of 80 minutes.

Of course also higher concentrations of both penicillin G and sulphadiazine can be detected.

5

Table 2: Sensitivity of test type C at different spore concentrations.

	Test duration (minutes)	140	123	105	90	80
10	spores/ml medium	1.3×10^7	2.5×10^7	5.1×10^7	1.0×10^8	2.0×10^8
	milk without antibiotics	-	-	-	-	-
	milk with 3 ppb penicillin G	+	+	+	+	+
15	milk with 100 ppb sulphadiazine	+	+	+	+	-

EXAMPLE VIII

20

Detection of antibiotics in milk (test type D)

In this experiment *Streptococcus thermophilus* T101 (DSM 4022) was used for the determination of antibiotics in milk.

Streptococcus thermophilus T101 was grown using methods known per se.

25

An amount of a solution of bromocresolpurple to give a final concentration of 45 mg per 1000 ml was added to fresh cow's milk.

To part of the fresh cow's milk also an amount of penicillin G of 4 ppb was added.

30

The control sample contained no penicillin G.

- 18 -

To the milk samples different amounts of *Streptococcus thermophilus* T101 were added to give final concentrations between 10^7 and 10^9 CFU per ml.

Immediately after preparation of the samples sterile tubes with a diameter of about 9 mm were filled with 0.3 ml of each milk sample (duplo).

5 The tubes were immediately placed in a waterbath kept at 42°C.

Observations were made after 30 minutes to 2 hours and 30 minutes.

The time at which the colour of the agar medium became yellow was determined. As is shown in table 3, the test duration correlated with the amount of cells added to the milk without influencing the sensitivity of the test in a negative way.

Even at a test duration of 80 minutes penicillin G can be detected at a concentration of 4 ppb or higher.

Table 3: Sensitivity of test type D at different cell concentrations.

15

20

25

Test duration (minutes)	145	125	103	80
cells/ml medium	3.0×10^7	6.0×10^7	$1,3 \times 10^8$	$2,5 \times 10^8$
milk without antibiotica	-	-	-	-
milk with 4ppb penicillin G	+	+	+	+

Claims

1. A method for detecting one or more antibacterial compounds in a liquid sample which comprises the steps of:

5 (i) contacting said sample with a test microorganism selected from thermophilic strains of *Bacillus* and *Streptococcus*, prior to addition of said sample growth of said test organism prevented;

(ii) providing conditions for said test microorganism whereby in the absence of antibacterial compounds growth of said test microorganism can
10 occur and can be detected; and

(iii) detecting the growth of said test microorganism after a time sufficient for detection of test microorganism growth in the absence of antibacterial compounds, but during which period detectable inhibition of test microorganism growth is maintained by a known amount of an antibacterial
15 compound in a test sample;
whereby said test organism is present at above 10^7 colony-forming units (CFU)/ml.

2. A method according to claim 1 wherein the test microorganism is
20 dispersed in a solidified medium preferably in a solidified agar medium.

3. A method according to claim 1 wherein the test microorganism prior to addition of the test sample, is in concentrated form preferably in dried form.

25 4. A method according to any one of claims 1-3 suitable for detecting in less than 2.5 hours Penicillin G at 0.003 I.U. per ml.

5. A method according to any one of claims 1-4 wherein prior to addition of the sample, said test microorganism is present in said agar medium
30 at 2×10^7 to 3×10^8 CFU/ml, or in case the test microorganism is not present in

- 20 -

said medium, the test microorganism is after the addition of the sample present at 2×10^7 to 3×10^8 CFU/ml sample.

5 6. A method according to claim 1 whereby the growth of said test microorganism is detected by colour change of the agar medium resulting from presence of an acid-base or redox indicator.

10 7. A method according to any one of claims 1 to 6 wherein the test microorganism employed is a strain of *Bacillus Stearothermophilus*, preferably wherein *Bacillus Stearothermophilus* var. *calidolactis* C953 is employed, or wherein the test organism employed is a strain of *Streptococcus thermophilus*, preferably wherein *Streptococcus thermophilus* T101 is employed.

15 8. A method according to claim 2 wherein said agar medium has been formed from an agar solution containing sufficient nutrients to support growth of said test organism.

20 9. A method according to claim 2 wherein the nutrient source for said test organism, or part thereof, is applied to the solidified agar medium.

10. A method according to claim 2 wherein said agar medium contains at least one substance to change the sensitivity of the test organism to antibacterial compounds.

25 11. A method according to claim 2 wherein 10 μ l - 3 ml agar medium is used.

30 12. A test unit for carrying out a method according to claim 2 comprising a container holding a solidified agar medium having dispersed therein a test organism selected from thermophilic strains of *Bacillus* and *Streptococcus* at above 10^7 CFU/ml, said agar medium optionally additionally

- 21 -

containing nutrients for growth of said test organism and/or one or more components for detection of said growth.

13. A test unit as claimed in claim 12 wherein said container is a tube or
5 a block with tubular indentations.

Figure 1
Test duration versus spore-concentration


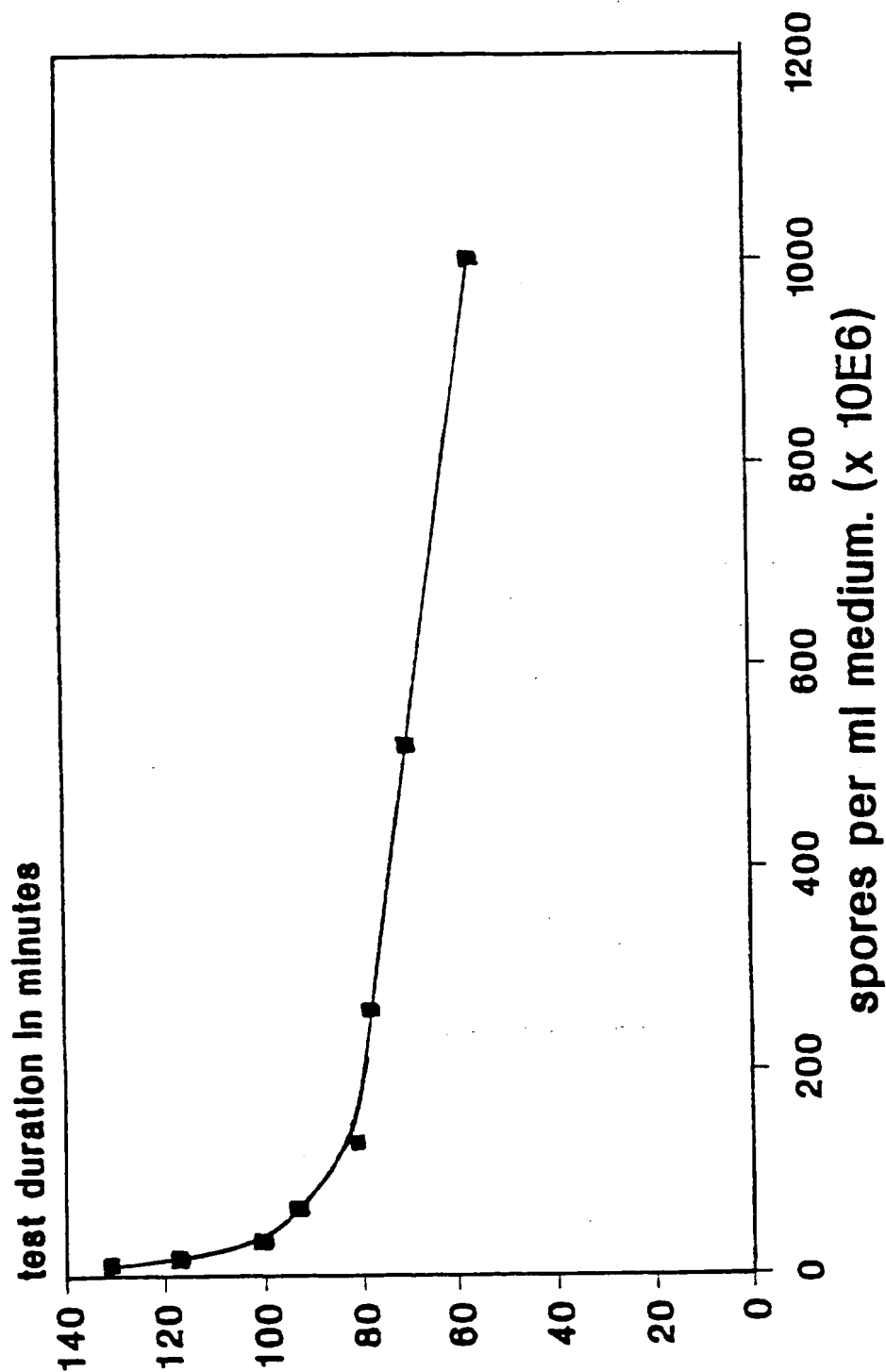


Figure 2
Test duration versus spore-concentration

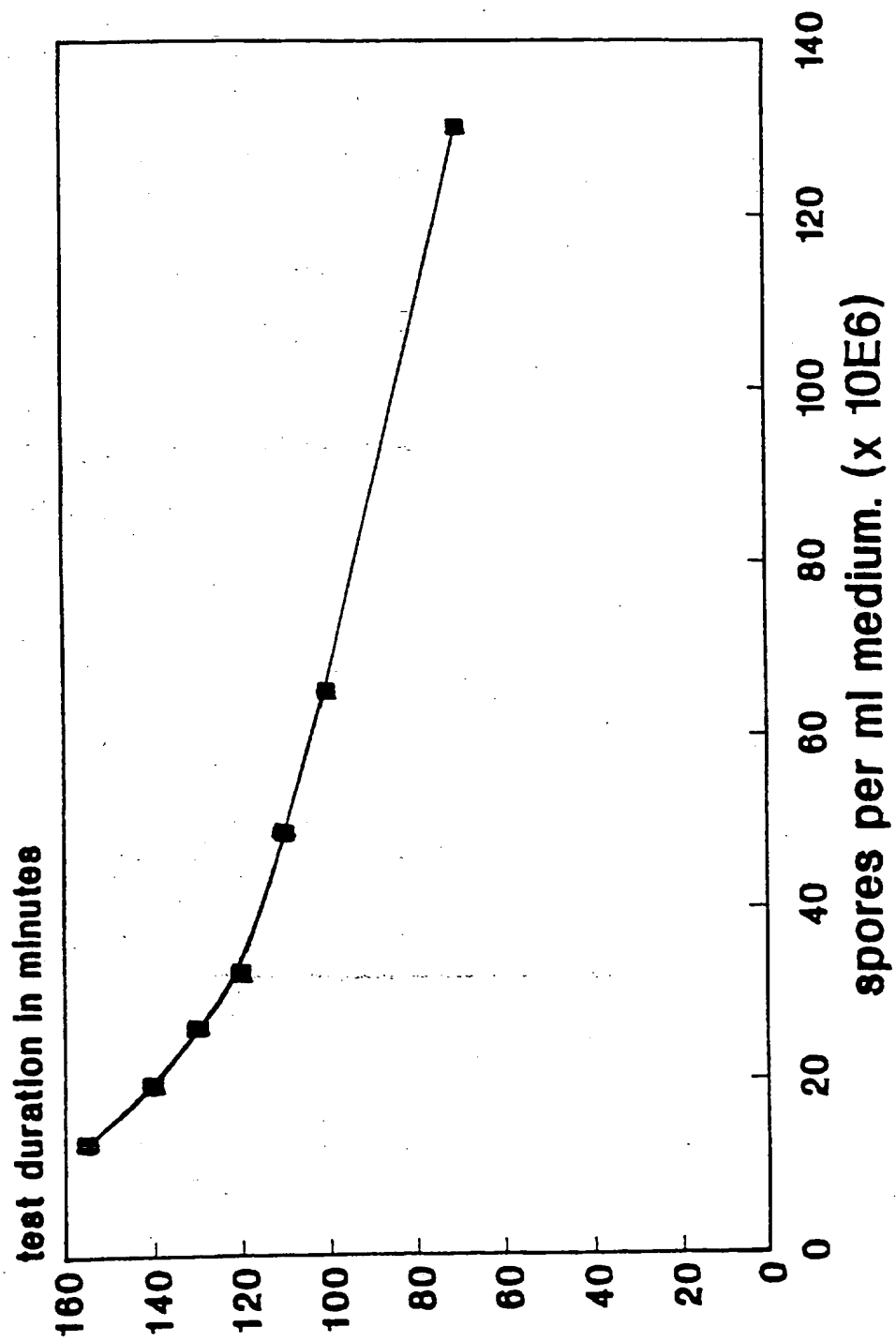

INTERNATIONAL SEARCH REPORT

International Application No
PC, EP 96/00488A. CLASSIFICATION OF SUBJECT MATTER
IPC 6 C12Q1/18 C12Q1/04

According to International Patent Classification (IPC) or to both national classification and IPC

B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)

IPC 6 C12Q

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practical, search terms used)

C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	US,A,4 946 777 (S.A. LAMÉRUS ET AL.) 7 August 1990 cited in the application see column 4, line 5 - line 22; example 1 AND 2 -----	1-13

☐ Further documents are listed in the continuation of box C.☒ Patent family members are listed in annex.

* Special categories of cited documents:

- *A* document defining the general state of the art which is not considered to be of particular relevance
- *E* earlier document but published on or after the international filing date
- *L* document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)
- *O* document referring to an oral disclosure, use, exhibition or other means
- *P* document published prior to the international filing date but later than the priority date claimed

T later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention

X document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone

Y document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.

A document member of the same patent family

Date of the actual completion of the international search

29 April 1996

Date of mailing of the international search report

30.05.96

Name and mailing address of the ISA

European Patent Office, P.B. 5818 Patentlaan 2
NL - 2280 HV Rijswijk
Tel. (+31-70) 340-2040, Tx. 31 651 epo nl,
Fax (+31-70) 340-3016

Authorized officer

Van Bohemen, C

INTERNATIONAL SEARCH REPORT

Information on patent family members

International Application No

PCT, EP 96/00488

Patent document cited in search report	Publication date	Patent family member(s)	Publication date
US-A-4946777	07-08-90	GB-A- 1467439	16-03-77
		AR-A- 210724	15-09-77
		AU-B- 6953574	04-12-75
		BE-A- 815759	02-12-74
		CA-A- 1040517	17-10-78
		CH-A- 611033	15-05-79
		DE-A- 2425999	19-12-74
		FR-A,B 2231966	27-12-74
		JP-C- 1206835	11-05-84
		JP-A- 50033897	01-04-75
		JP-B- 58032879	15-07-83
		NL-A- 7407251	03-12-74
		SE-B- 418621	15-06-81
		SE-A- 7407176	02-12-74
		US-A- 3941658	02-03-76

THIS PAGE BLANK (USPTO)

**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ **BLACK BORDERS**
- ☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**
- ☒ **FADED TEXT OR DRAWING**
- ☒ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**
- ☐ **SKEWED/SLANTED IMAGES**
- ☐ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**
- ☐ **GRAY SCALE DOCUMENTS**
- ☐ **LINES OR MARKS ON ORIGINAL DOCUMENT**
- ☐ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**
- ☐ **OTHER:** _____

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.

THIS PAGE BLANK (USPTO)